

HOW TO INSTALL YOUR LARES REPLACEMENT TURBINE CARTRIDGE STEP BY STEP.

If replacements cannot be made successfully, please return the entire handpiece to Lares Research for factory repair.




Should you need information about the operation, maintenance, or service of your Lares product, contact us directly. Lares Research is here to help you.

Within the U.S.A. 1-888-333-8440, Ext. 1 Worldwide 530-345-1767, Ext. 1

Outside the USA contact your Authorized Lares Distributor.

PROCEDURE

1. Wear exam gloves.
2. Secure an expended standard length bur in chuck.
3. Mark head cap and handpiece with a permanent marker before opening. This is critical to assure head cap is tightened properly when replacement turbine is installed (Figure 1).
4. Engage head cap with head cap wrench as shown below. (Figure 2A or 2B, depending on model). Grasp cap wrench firmly. Turn counterclockwise with inward pressure until threads are disengaged. Note: The same force required to loosen cap should be used when tightening cap during reassembly.
5. Making sure the cap is fully disengaged, press bur down against flat surface to remove old cartridge and cap from head (Figure 3). Remove head cap assembly from turbine cartridge. Discard worn turbine cartridge and retain head cap. Locate and retain both wavesprings, one between the cap and turbine bearing and the other one between the head bore and turbine bearing (below o-rings) (Figure 7). Note: If ball bearing is destroyed we recommend that the outer ring of the ball bearing, which remained in the turbine housing or in the cap, is removed by an authorized Lares repair facility by means of a special tool.
6. Flush interior of head and head cap assembly with alcohol and clean with cotton swab. Flush again and inspect for cleanliness. Repeat if necessary. Note: Please do not open the new turbine pouch until the handpiece has been prepared as above. This will assure maximum life and performance from the turbine. The pouch and turbine have been manufactured and sealed within a clean zone, which protects the precision bearings from possible contamination.
7. Apply Lares conditioner to the entire area of the head and head cap interior. Verify position of wave springs, one in head cap and one in head bore (beyond o-ring) (Figure 7). Remove the new turbine from the pouch.

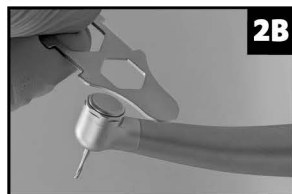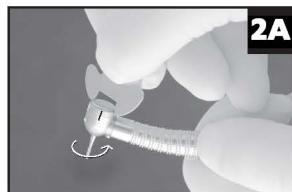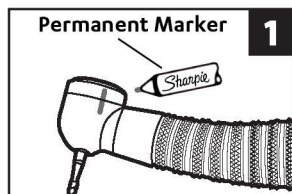




ProStyle SF/Legacy 5

Xtreme Power (XP)

PROCEDURE CONTINUED

8. Insert turbine cartridge into head cap assembly until fully seated. Then carefully insert entire assembly (bur end first) into turbine head (Figure 4) until cap rests against back of turbine head. Turn counterclockwise with light pressure until seated as indicated by "click" sound (Figure 5). Then turn clockwise with light inward pressure until fully engaged in head.
 9. Engage head cap with head cap wrench as shown in Figure 2A or 2B, depending on model. With inward pressure, turn head cap wrench clockwise until fully tightened. Use same force to tighten as was used to loosen it. Tighten the head cap with the wrench until it is tight and the marked lines meet.
 10. Double check to make sure that the head cap is fully tightened so that the mark on the head cap and head line up.
-  **WARNING:** If cap is tightened improperly (too loose) it may result in the turbine ejecting itself during use.
11. With a bur in place connect and run the handpiece at half speed for a minimum of one minute. Wipe off any visible lubricant.

Your handpiece is now ready to operate.

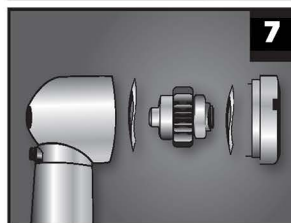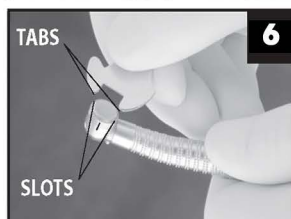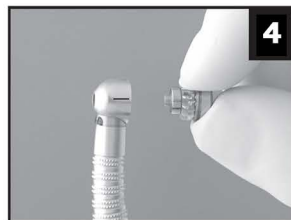

LIMITED WARRANTY

Lares replacement turbine cartridges with push button bur changing are warranted against defects in materials and workmanship for a period of one year from date of shipment.

1. The handpiece must be operated and maintained in accordance with procedures outlined in the instructions.
2. The handpiece must not be subjected to abuse or neglect.
3. The handpiece must not have been repaired or disassembled by anyone other than Lares Research or your Authorized Lares Distributor.

Lares Research will repair or replace, at its discretion, without charge any defective parts covered by this warranty provided the handpiece is returned to the factory. (Outside the U.S.A. return to your Authorized Lares Distributor). Lares Research makes no other warranties, expressed or implied.

lares[®]
RESEARCH
Precision Performance

295 Lockheed Avenue
Chico, CA 95973 U.S.A
530-345-1767 • Fax 530-345-1870
www.laresdental.com