


INTRODUCING THE WORLD'S FIRST TURBINE INTEGRATED CARIES DETECTION

lares[®]
RESEARCH
60 Years of Quality



FLUORESCENCE™ HD

- 405 nm light emitted from turbine or contra angle causes caries to fluoresce red/orange and healthy tooth light green.
- Dentist removes only orange/red diseased tissue, minimizing removal of healthy tooth.
- F.A.C.E. (Fluorescence Aided Caries Excavation) - is more accurate for detecting remaining decay than probing or caries indicator dye, minimizing risk of reinfection.
- Eliminating the need for caries indicator dyes or probing saves chair time.

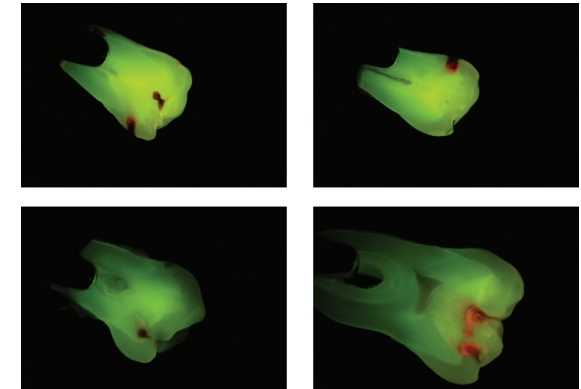

Fluoresce HD Provides Constant Visual Differentiation Of Remaining Caries Vs. Healthy Tissue.



Fluoresce HD can be used for decay removal with both Lares Fiber Optic Turbines and Fiber Optic Lowspeeds



LARES BOOTH:

U.S. Pavillion, Hall 4.2, Aisle H-069
Contact Christian Godoy: cgodoy@laresdental.com
Mobile: 1-530-717-3145 www.laresdental.com

INTRODUCING DENTISTRY'S FIRST 5 YEAR TURBINE WARRANTY

lares[®]
RESEARCH
60 Years of Quality

Eliminate Your Repair Costs For 5 Years With New Lares LEGACY 5 Turbines.

Now you can eliminate the expense and hassle of handpiece repair charges for five years with the new Lares Legacy 5 Turbines. Say good-bye to outrageous repair charges, and hello to exceptional turbine performance, legendary Lares service and never-before convenience!



LEGACY 5 Highspeeds have the Features and Performance you Need

- Dentistry's only 5 year turbine warranty!
- Choice of Standard or Mini head sizes.
- Exceptional cutting power (standard: 20 W, mini: 18 W)
- DynamicPrecision™ balancing for whisper quiet operation (Standard head 59.9 dBA)
- SteelGrip™ push button chucking system for long lasting bur retention and cutting precision.
- Solid rod cellular optics for bright, long lasting illumination (optional).
- KaVo* MULTiflex* compatible 360° swivel and quick disconnect.
- Stainless steel construction with non-slip PVD coated grip.
- Long life precision ceramic bearings.



*Legacy 5 Mini Head
For Excellent Visibility and
Accessibility*

LARES BOOTH:

U.S. Pavillion, Hall 4.2, Aisle H-069

Contact Christian Godoy: cgodoy@laresdental.com

Mobile: 1-530-717-3145 www.laresdental.com



*Legacy 5 Standard Head
For Crown and Bridge
With 4 Port spray*

*Kavo and MULTiflex are trademarks of Kavo Dental GMBH.