


# HANDPIECE DEPLOYED CARIES DETECTION




**NEW!**

## FLUORESCENCE™ HD

- During caries removal, the 405 nm (UV) LED light causes caries to fluoresce orange/red and healthy tooth light green.
- Dentist removes only orange/red caries, minimizing the removal of healthy tooth (minimally invasive).
- When all orange/red decayed tissue is removed, preparation is complete and no need to change instruments which substantially reduces chair time.
- F.A.C.E. (Fluorescence Aided Caries Excavation) is more accurate at identifying remaining decay than alternative methods, minimizing the risk of reinfection.
- Fluoresce also dramatically accentuates visual differences between tooth and restorative materials (i.e. resins, cements, etc.) further minimizing the removal of healthy tooth and saving chair time.

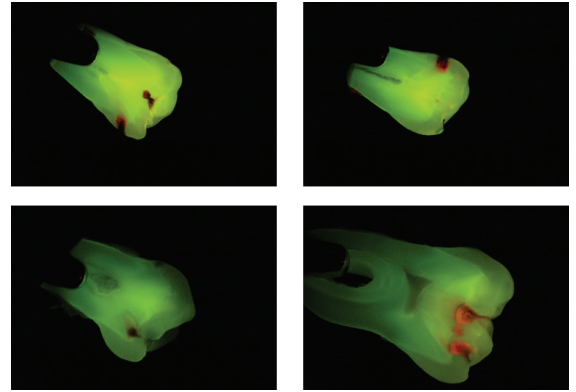

*Fluoresce HD Provides Clear Visual Differentiation Of Remaining Caries Vs. Healthy Tissue.*

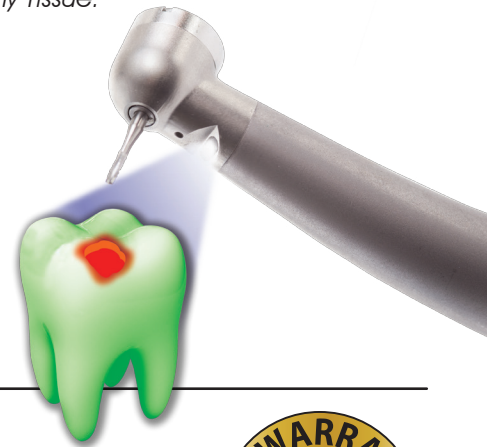

**NEW!**

## LEGACY 5 TURBINES: BREAKTHROUGH 5 YEAR WARRANTY



### LEGACY 5 Highspeeds have the Features and Performance you Need

- Exceptional cutting power (standard: 20 W, mini: 18 W)
- DynamicPrecision™ balancing for whisper quiet operation (standard head 59.9 dBA)
- SteelGrip™ push button chucking system for long lasting bur retention and cutting precision.
- Solid rod cellular optics for bright, long lasting illumination (optional).
- KaVo MULTiflex\* compatible 360° swivel and quick disconnect.
- Stainless steel construction with non-slip TactiSense™ grip.
- Long life precision ceramic bearings.



**Legacy 5**  
Standard  
Head with  
4 Port Spray



**Legacy 5 Mini Head**  
For Excellent Visibility and  
Accessibility

## SEEKING INTERNATIONAL DISTRIBUTORS

Contact Christian Godoy: [cgodoy@laresdental.com](mailto:cgodoy@laresdental.com)  
Mobile: 1-530-717-3145 [www.laresdental.com](http://www.laresdental.com)

\*KaVo and MULTiflex are trademarks of KaVo Dental GmbH